

FIGURE 9.24 As described in the legend to Fig. 9.22 but for m = 100. Note that the distribution is centered at a mean time $t \approx 3.56$ s, where $t = (m-1)/r \approx 100/r$.

100 events is much more "stable" (relative to its mean value) than between every second event. As can be seen from Eq. (9.24) the distributions for m > 1 have a maximum $(dq_m/dt = 0)$ at

$$t = \frac{m-1}{r}. (9.27)$$

Thus, from the location of the peak in the distribution we can obtain the average rate. We find that for the data shown in Figs. 9.23, and 9.24

$$m = 3$$
, $t_{\text{max}} = 0.073 \text{ s}$, $r = 27.5 / \text{s}$
 $m = 100$, $= 3.56 \text{ s}$, $= 27.7 / \text{s}$.

Furthermore, a fit to the exponential for m = 1 (see Fig. 9.22) yields $t_{1/e} = 3.50 \times 10^{-2}$ s, or r = 28.6/s in agreement with the average rate.

9.4.3. Measurement of the Mean Life of the Muon

The muon is not stable but decays into an electron, a neutrino, and an antineutrino:

$$\mu^{+} \rightarrow e^{+} + \nu_{e} + \overline{\nu}_{\mu}$$

$$\mu^{-} \rightarrow e^{-} + \overline{\nu}_{e} + \nu_{\mu}.$$
(9.28)

The mean life, or lifetime, (i.e., the inverse of the decay rate) for this process is of order 2.2 µs, and thus the decay is easily detectable for muons at rest. The neutrinos are not observable but the electron (or positron)²⁴ is energetic enough to give a clear signal of the decay. The mass of the muon is

$$m_{\mu} = 105.65 \,\mathrm{MeV/c^2},$$

approximately 200 times the electron mass. The maximum energy of the electron occurs when the two neutrinos recoil against it as shown in Fig. 9.25a. This corresponds, in the rest frame of the muon, to

$$E_e(\text{max}) \simeq \frac{1}{2} m_\mu c^2 = 53 \text{ MeV}.$$

The energy spectrum of the electrons from muon decay is shown in Fig. 9.25b.

The long lifetime for muon decay indicates that the decay does not proceed through the strong (nuclear) interaction but rather through the weak interaction responsible for the " β -decay" of nuclei. However, the process of Eq. (9.28) is very important because it involves only leptons (no strongly interacting particles participate) and thus can be used unambiguously to

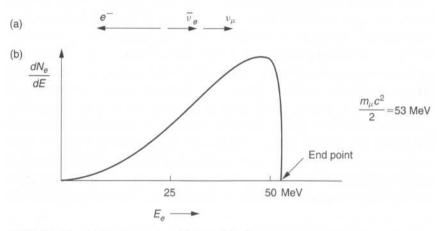


FIGURE 9.25 (a) Configuration of the particles in μ -decay for obtaining the maximum electron energy. (b) The energy spectrum of the electrons from μ -decay.

²⁴To save words we will speak only of the electron even though we mean either e^- or e^+ .

calculate the Fermi weak interaction constant G_F . The mean life of the muon is given by

$$\frac{1}{\tau} = \frac{1}{\hbar} \frac{G_{\rm F}^2}{(\hbar c)^6} \frac{(m_\mu c^2)^5}{192\pi^3}.$$
 (9.29)

Precise measurements of the muon mean life yield

$$\tau_{\mu} = (2.19703 \pm 0.00004) \times 10^{-6} \text{ s}$$
 (9.30a)

and through Eq. (9.29) the value of the Fermi constant²⁵

$$\frac{G_{\rm F}}{(\hbar c)^3} = 1.1664 \times 10^{-5} \,\text{GeV}^{-2}.$$
 (9.30b)

We will measure the decay of muons that have come to rest in the liquid scintillator tank. Muons lose approximately 2 MeV of energy for each gram per squared centimeter of material that they traverse. Thus we expect that muons entering the 35-cm-high liquid scintillator tank with energy $E_{\mu} \lesssim 50$ MeV will stop in the tank. The fraction of muons that do stop is of order 0.3% of the flux going through the tank. Thus, the stopping rate is $R_{\rm S} = 0.077/{\rm s}$, or 4.6 muons/min. This is adequate to obtain good statistical accuracy for the mean life value in a reasonable time interval.

The experimental arrangement is shown in Fig. 9.26. When a muon enters the tank the PMT gives a signal, which is amplified and then discriminated. This pulse is used to start a "time-to-amplitude converter" (TAC) circuit. If the muon stops in the tank, then the decay electron will give a second signal within a time interval of a few mean lives. The second pulse is used to stop the TAC, and the time interval between start and stop is directly read out. The 60-ns delay in the start signal is to make sure that no pulse will be on the stop line when the start arrives. Commercial electronic modules can be used to achieve the logic indicated in Fig. 9.26. A GT200 computer card, designed by Professor D. Hartill of Cornell University, 26 performed the TAC functions and stored the data in a file in the computer memory. If no stop arrives within $\Delta t = 25 \,\mu s$, which corresponds to ~ 10

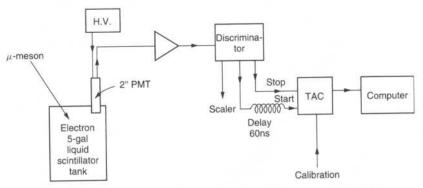


FIGURE 9.26 Block diagram of the electronics for measuring muon decay.

mean lives, the TAC is reset and the start pulse ignored. To calibrate the TAC one applies a fixed frequency (oscillator) signal to the discriminator input.

If the singles rate is too high, then the stop pulse may not be due to the decay of the muon that started the TAC, but to a different muon entering the tank. We call such events "accidental stops," and we can estimate their rates as follows. The singles rate is r=25/s, so that using the Poisson distribution of Eq. (9.23) for n=2 and $t=\Delta t=25$ μs we find for the accidental rate

$$R_{\rm a} = \frac{P_{\rm a}(n=2, \Delta t)}{\Delta t} = 7.8 \times 10^{-3} \,\mathrm{s}^{-1}.$$
 (9.31)

This is ten times smaller than the stopping rate R_S , and does not affect the determination of the mean life as discussed later.

Data obtained by a student are shown in Fig. 9.27. The data were accumulated over five days and yielded $N_{\rm S}=32,000$ stops in 6921 min. The very early events, $t<0.25~\mu \rm s$, were discarded, leaving a sample of 30,069 events displayed in 100 bins each 0.25 $\mu \rm s$ wide. The data for $t\lesssim 5~\mu \rm s$ show an exponential drop-off, as expected, and in this region are well fitted by

$$N(t) = N_0 e^{-t/\tau}$$
 0.25 < t < 5 \mus.

In contrast the data for late times, $t>15~\mu s$, are flat and are well fitted by a constant

$$N(t) = C$$
 5 < t < 25 \mus.

²⁵Note that in contrast to the fine structure constant $\alpha = e^2/\hbar c$, which is dimensionless, $G_F/(\hbar c)^3$ has dimensions of inverse energy squared. In fact $G_F/(\hbar c)^3$ has the approximate value of $1/(M_W c^2)^2$, where M_W is the mass of the vector bosons that mediate the weak interactions.

²⁶One can now purchase commercial versions of TAC cards to perform the required functions.

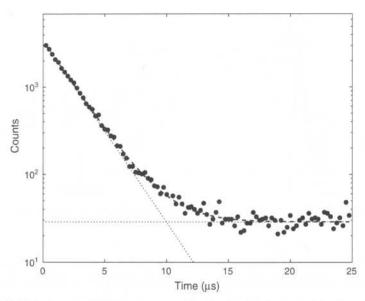


FIGURE 9.27 Data for 30,000 muon stops. The bin size is $0.25\,\mu s$, and the fit to the data including an exponential decay and a constant background term are shown.

A combined fit²⁷ of the form

$$N(t) = N_0 e^{-t/\tau} + C (9.32)$$

yields $\tau = 2.088 \pm 0.016 \,\mu\text{s}$, $N_0 = 3410$, and C = 28.8, and an excellent $\chi^2 = 0.909$ per degree of freedom. The contributions of the two terms of the fit are also indicated by the dashed lines in the graph.

We briefly discuss the background level. Since there are 100 channels, the total accidental count is $N_{\rm a}=2880$, and thus the accidental rate is $R_{\rm a}=N_{\rm a}/6921$ min = $6.9\times10^{-3}/{\rm s}$ in agreement with our estimate of Eq. (9.31). One recognizes that the background does not affect the measurement until

$$N_0e^{-t/\tau}\sim C.$$

This occurs when $t/\tau \sim 4.71$, which allows for a fairly long "lever arm" to determine τ_u .

Our value for the mean life is in close agreement with the accepted value as given in Eq. (9.30a). The agreement is even closer because the

measured value for τ_{μ} must be corrected for the following effect. When negative muons stop in matter, there is a finite probability that the μ^- will be absorbed by a proton in the nucleus, leading to a "capture" reaction:

$$\mu^- + Z \to (Z - 1)^* + \nu_{\mu}$$
.

Thus the effective mean life is shortened and given by

$$\frac{1}{\tau_e} = \frac{1}{\tau_u} + \frac{1}{\tau_c},$$

where $1/\tau_{\mu}$ and $1/\tau_{c}$ are the rates for decay and capture, respectively. As a result the observed mean life is shorter; for mineral oil (the capture occurs mainly on carbon nuclei) and for the μ^{-}/μ^{+} composition of cosmic rays this correction is approximately 4%. Therefore, the corrected measured value in this experiment is

$$\tau_{\mu} = 2.172 \pm 0.017 \,\mu \text{s}.$$
 (9.33)

The error shown in Eq. (9.33) is only statistical and does not include systematic effects, in particular any uncertainty in the TAC calibration.

9.5. γ-γ ANGULAR CORRELATION MEASUREMENTS

9.5.1. General Considerations

We will now discuss the measurement of the correlation in angle between two gamma rays that were emitted simultaneously from the same source. The origin of these gamma rays is frequently the cascaded decay of a nucleus, as in the case of ⁶⁰Ni (⁶⁰Co) already discussed in Chapter 8. (See Fig. 8.20.) We reproduce in Fig. 9.28 the decay scheme of this nucleus and note that the 1.333-MeV gamma ray follows the 1.172-MeV gamma ray, the lifetime of the intermediate state being only about 10^{-12} s, so that for all practical purposes the two gamma rays are coincident.

The fact that these two gamma rays are correlated in angle can be understood from the following general argument: the first gamma ray will have an angular distribution with respect to the spin axis of the nucleus; thus its observation at a fixed angle $\theta=0$ conveys information about the probability of finding the spin at some angle ψ with respect to the

²⁷See Section 8.6.2.